
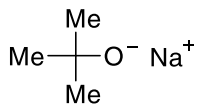
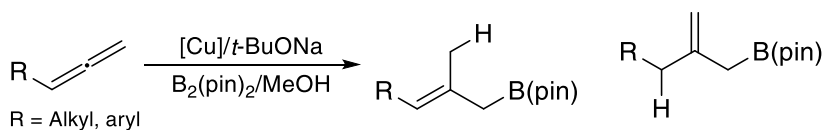
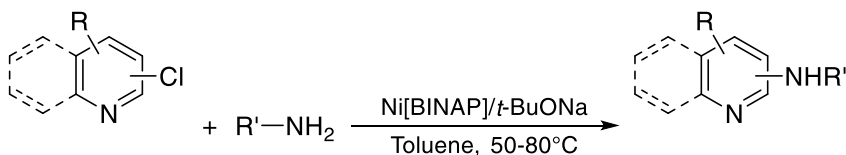
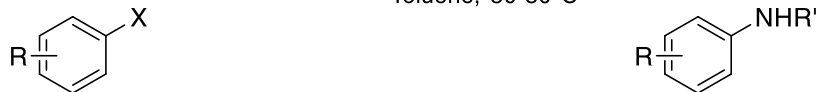


Catalog # 93-1022 CALLERY™ Sodium *tert*-butoxide, min. 98%


Technical Notes:

1. Base used in highly selective Cu-catalyzed hydroboration of allenes and 1,3-dienes.
2. Used in the Ni[BINAP]-catalyzed amination of aryl and heteroaryl chlorides and bromides with primary aliphatic amines.
3. Base used in the Cu-catalyzed cross-coupling of boronic esters with aryl iodides and carboration processes of alkynes and allenes.
4. Used in the enantioselective Cu-catalyzed 1,6-conjugate additions of propargyl and allyl groups.
5. Catalyst for transition-metal-free and redox-neutral lactamization of sp^2 C–H bonds with CO_2 .
6. Used in the base-promoted synthesis of *N*-substituted 1,2,3-triazoles via enaminone-azide cycloaddition involving Regitz diazo transfer.
7. Used as a pre-catalyst activator in Fe- and Co-catalysed hydroboration and hydrosilylation of alkenes and alkynes.
8. Used in ligand-controlled preparation of alkyl aryl ethers via Cu-catalyzed alkoxylation of (hetero)aryl halides.
9. Used in Ni/NHC-catalyzed asymmetric C–H alkylation of fluoroarenes with alkenes to generate enantioenriched fluorotetralins.
10. Used in Cu-catalyzed stereoselective coupling of terminal alkynes and α -bromo carbonyls to generate functionalized *E*-alkenes.
11. Used in Cu-catalyzed coupling reactions with (hetero)aryl chlorides and Bromides to generate α -(hetero)aryl nitriles.
12. Used in Pd-catalyzed enantioselective α -carbonylative arylation for facile construction of chiral spirocyclic β,β' -diketones.


Tech Note (1)
Ref. (1)Tech Note (2)
Ref. (2)Tech Note (3)
Ref. (3)

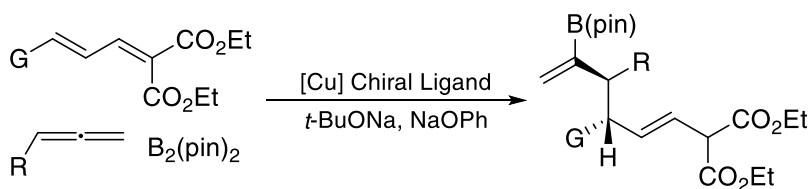


Tech Note (4)
Ref. (4)

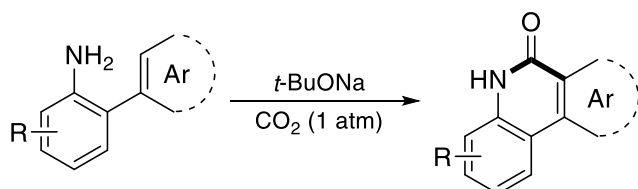


Tech Note (5)
Ref. (5)

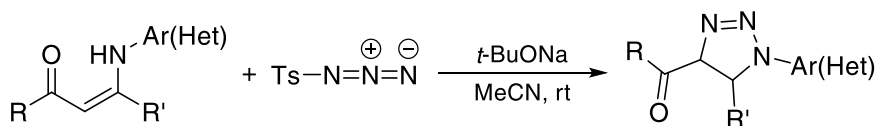


Tech Note (6)
Ref. (6)

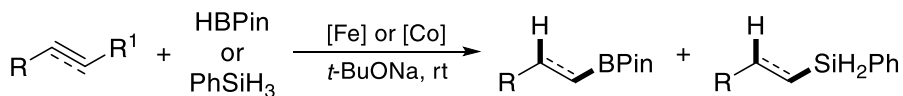

Tech Note (7)
Ref. (7)

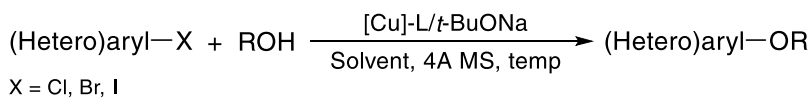


Tech Note (8)
Ref. (8)

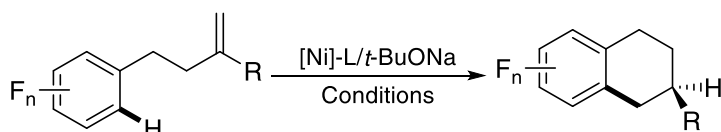


Tech Note (9)
Ref. (9)


$n = 2, 3, 4$; $R = \text{alkyl, aryl, alkenyl, or N}$

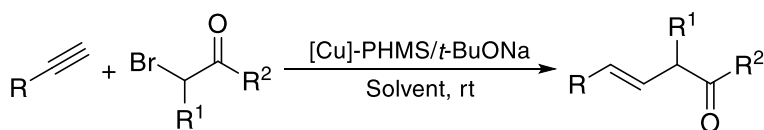

Tech Note (10)
Ref. (10)

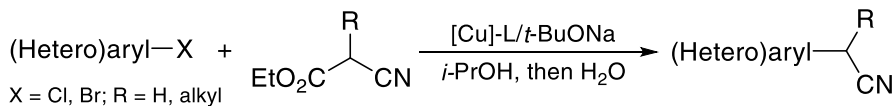

Tech Note (11)
Ref. (11)

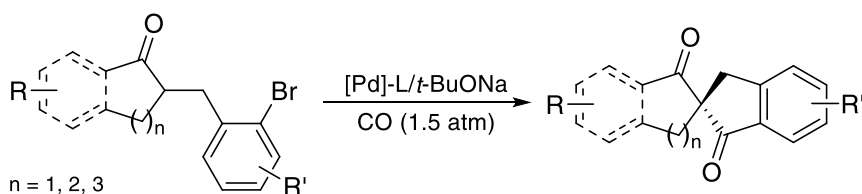

Tech Note (12)
Ref. (12)

References:

1. [Chem. Eur. J. 2013, 19, 7125.](#)
2. [J. Am. Chem. Soc. 2014, 136, 1617.](#)
3. [Angew. Chem. Int. Ed. 2014, 53, 3475.](#)

4. [Nature 2016, 537, 387.](#)
5. [Angew. Chem. Int. Ed. 2016, 55, 7068.](#)
6. [Org. Lett. 2016, 18, 6034.](#)
7. [Nat. Chem. 2017, 9, 595.](#)
8. [J. Am. Chem. Soc. 2019, 141, 3541.](#)
9. [Angew. Chem. Int. Ed. 2019, 58, 13433.](#)
10. [J. Am. Chem. Soc. 2021, 143, 7903.](#)
11. [Angew. Chem. Int. Ed. 2021, 60, 7082.](#)
12. [Angew. Chem. Int. Ed. 2021, 60, 9978.](#)